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FICHES TECHNIQUES - TECHNICAL SHEET

FICHES TECHNIQUES

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Fractal Plus	Ø 3.75 FPO	2	
	Ø 4.75 FPO	6	
Fractal Lift	Ø 3.75 OF1	10	
	Ø 4.75 OF1	14	
	Ø 3.75 MF1	18	
Fratex	Ø 3.30 H10 & H12 OFX	22	
	Ø 3.30 H16 OFX	26	
Fratex Oblik	Ø 3.30 MFX	30	
Pterygoïdien	Ø 4.75 MF4	34	

TECHNICAL SHEET

DOWNLOAD ALL 

Fractal Plus	Ø 3.75 FPO	4	
	Ø 4.75 FPO	8	
Fractal Lift	Ø 3.75 OF1	12	
	Ø 4.75 OF1	16	
	Ø 3.75 MF1	20	
Fratex	Ø 3.30 H10 & H12 OFX	24	
	Ø 3.30 H16 OFX	28	
Fratex Oblik	Ø 3.30 MFX	32	
Pterygoïdien	Ø 4.75 MF4	36	

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
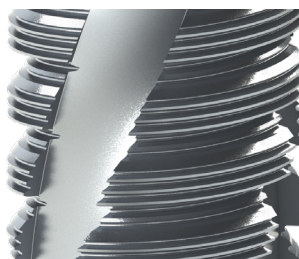
Roughness : Machined smooth area $0.6 \mu\text{m} < SA < 0.8 \mu\text{m}$


INDICATIONS

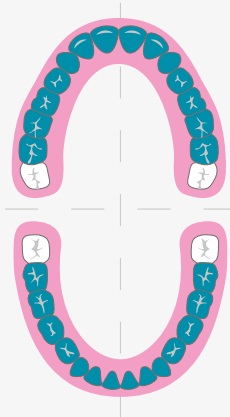- For all sectors replacement
- Prosthesis stabilization removable or fixed (unitary or interleaved). The volume and bone quality must be sufficient to support the implants.

SUPPLIED AS STANDARD WITH

- Internal octagon $\emptyset$ 4.50mm breakable cover screw

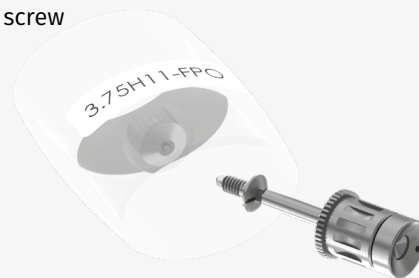Visy Implant

55, rue Uranus - Z.A.C. Altaïs - 74650 CHAVANOD (France) - Tel.: +33 (0) 450 450 498

@ support.technique@visyimplant.com victoryimplant.com / visyimplant.com


victoryimplant.com / visyimplant.com

MADE IN FRANCE

TECHNICAL SHEET Ref. 3.75-FPO

DRILLING PROTOCOL

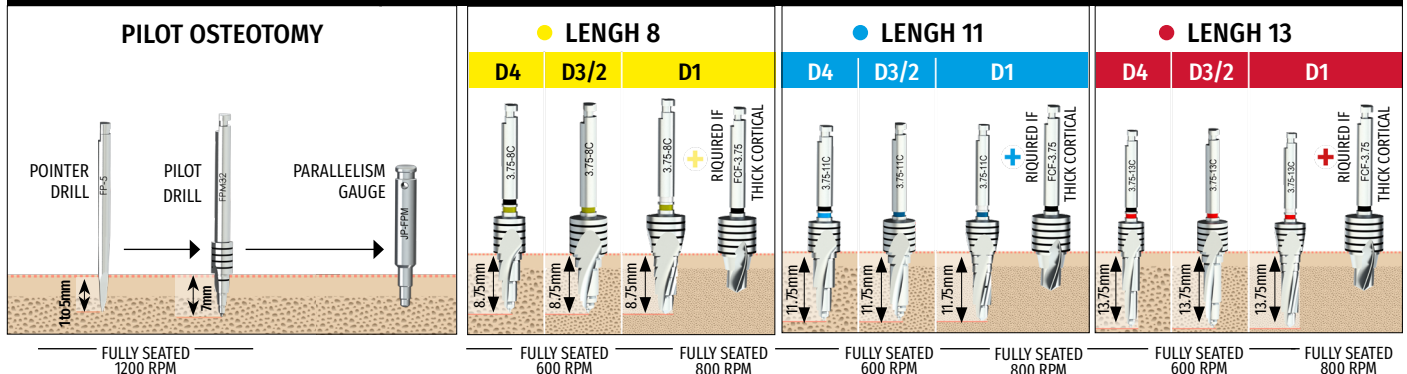
Drilling recommendation: Under continous irrigation.
Drilling torque recommendation: 5 to 10 N.cm. If the drill bit stalls at 10 N.cm, plan to replace it.

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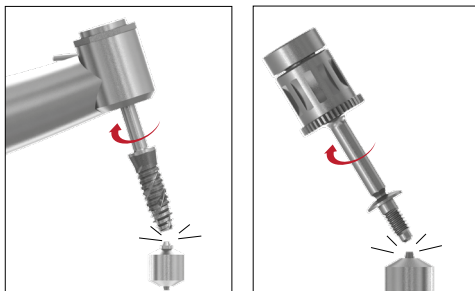
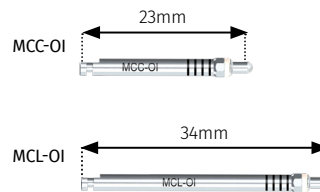
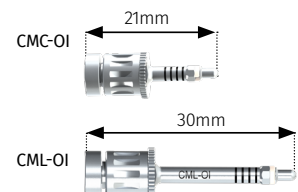
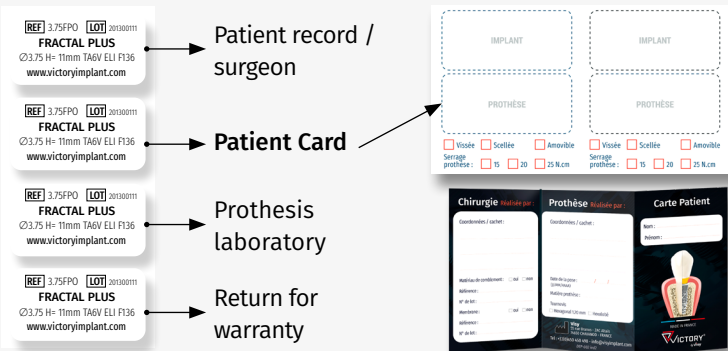
VICTORY WARRANTY

- Victory® products have a 10 years warranty, free from manufacturing defects.
- The buyer remains solely responsible for the use he makes of the products.
- Every compasant replacements are left to Visy discretion.

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
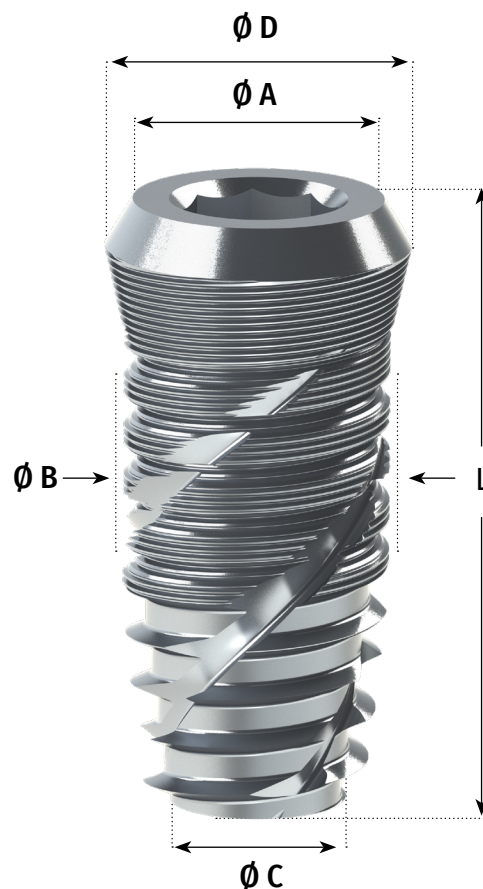
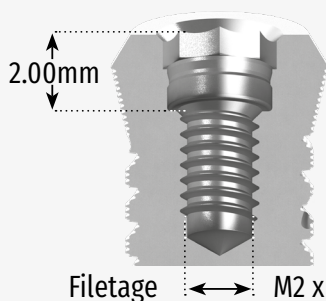
INDICATIONS

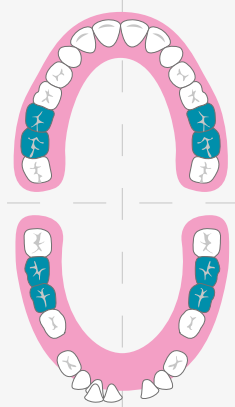- For the mollars replacement
- Prosthesis stabilization removable or fixed (unitary or interleaved). The volume and bone quality must be sufficient to support the implants.

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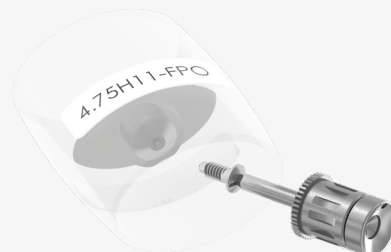
MADE IN FRANCE

TECHNICAL SHEET Ref. 4.75-FPO

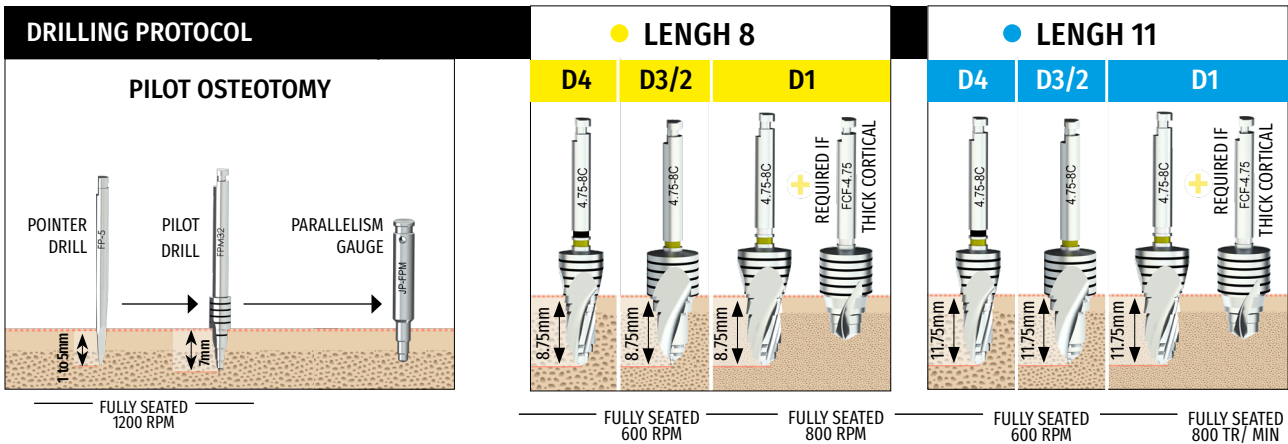i **Drilling recommendation:** Under continuous irrigation.
Drilling torque recommendation: 5 to 10 N.cm. If the drill bit stalls at 10 N.cm, plan to replace it.

! The drill tip has a length of 0.75 mm. This value must be taken into account during drilling in addition to the implant length. The values D1 to D4 indicate the bone densities for which the drill can be used. Do not use a rotational speed (in revolutions per minute) outside the indicated range.

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

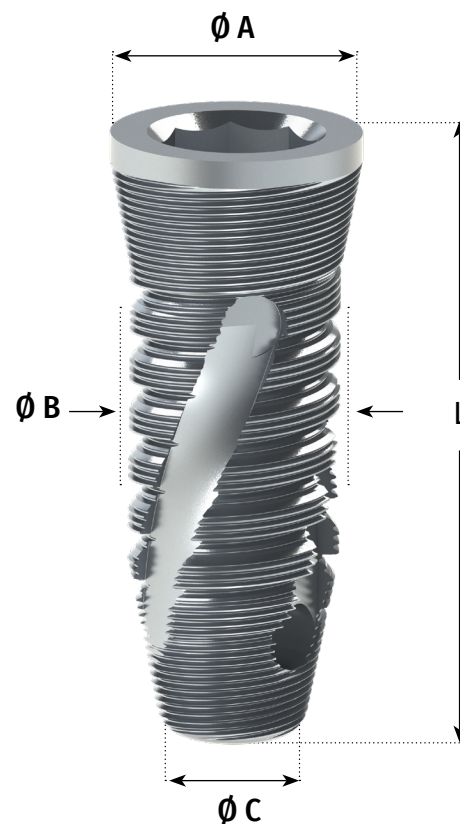
INTERNAL OCTAGON PLATFORM

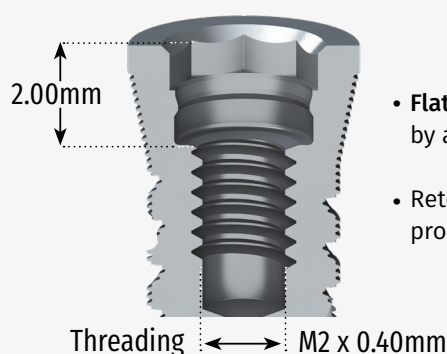- Flat indexed connection by a 2.6mm octagon
- Retention groove for centering prosthetic components

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
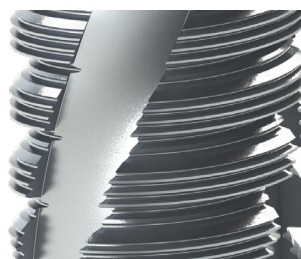
INDICATIONS

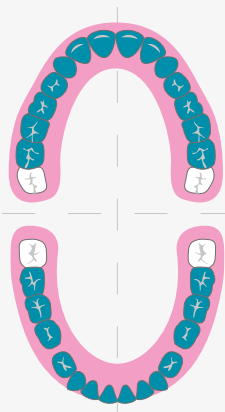- For the replacement in all sectors
- **Prosthesis stabilization** removable or fixed (unitary or interleaved). The volume and bone quality must be sufficient to support the implants.

SUPPLIED AS STANDARD WITH

- **Internal octagon $\varnothing$ 4.50mm** breakable cover screw

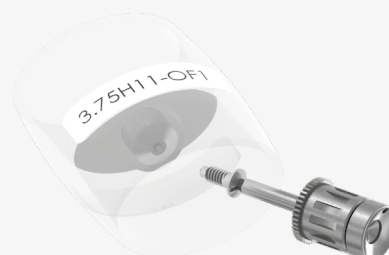Visy Implant

55, rue Uranus - Z.A.C. Altaïs - 74650 CHAVANOD (France) - Tel.: +33 (0) 450 450 498


@ support.technique@visyimplant.com  victoryimplant.com / visyimplant.com

MADE IN FRANCE

TECHNICAL SHEET Ref. 3.75-OF1

DRILLING PROTOCOL

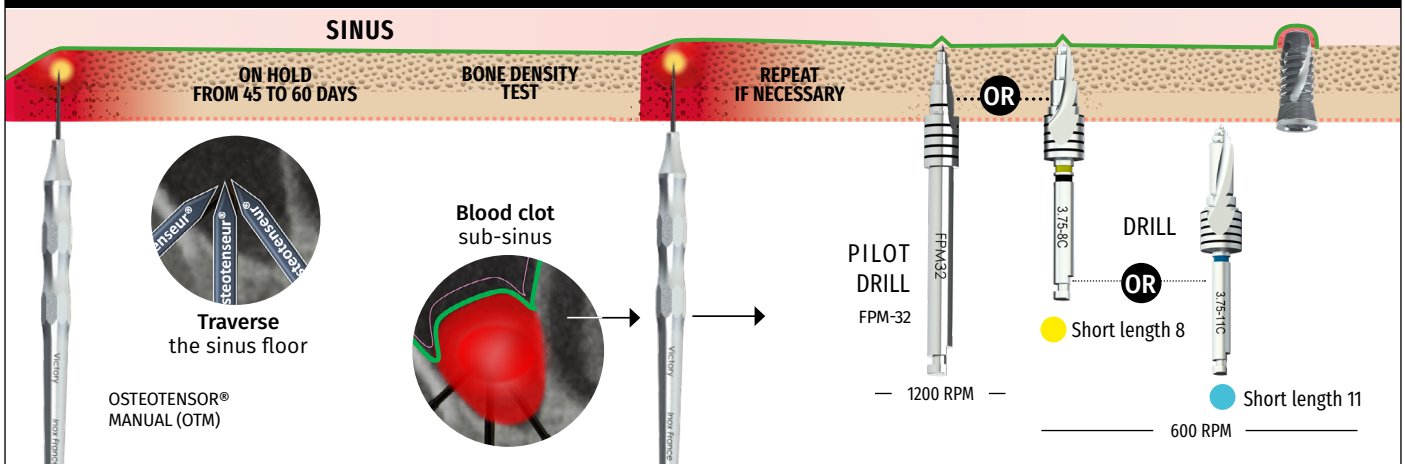Drilling recommendation: Under continuous irrigation.

Drilling torque recommendation: 5 to 10 N.cm. If the drill bit stalls at 10 N.cm, plan to replace it.

⚠ The drill tip has a length of 0.75 mm. This value must be taken into account during drilling in addition to the implant length.

The values D1 to D4 indicate the bone densities for which the drill can be used. Do not use a rotational speed (in revolutions per minute) outside the indicated range.

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


MONOBLOC PLATFORM

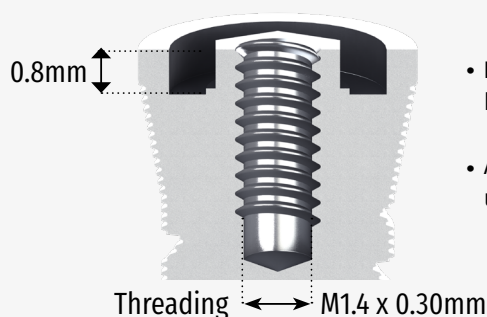- Indexed flat connection by a 2mm hexagon
- Allows axis correction up to 90° between two implants

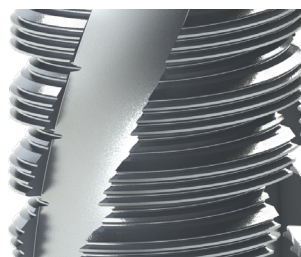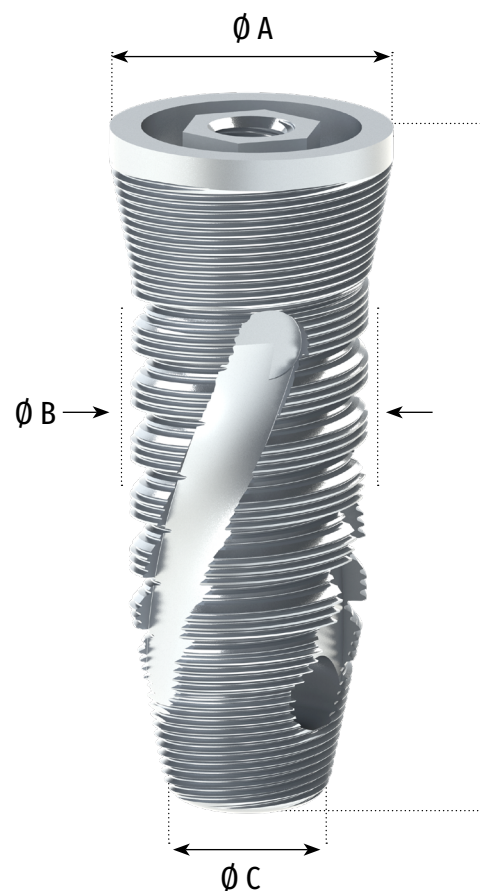Material


Titanium alloy TA6V ELI F136 - EN ISO 5832-3/ASTM F136

Surface state

VMSS (Victory Smooth Surface)

Roughness : Machined smooth area $0.6 \mu\text{m} < SA < 0.8 \mu\text{m}$

INDICATIONS

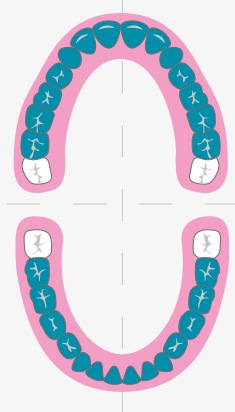- For all sectors replacement
- **Prosthesis stabilization** removable or fixed (unitary or interleaved). The volume and bone quality must be sufficient to support the implants.

SUPPLIED AS STANDARDS WITH

- **MONOBLOC $\varnothing$ 4.50mm** self breaking cover screw

Visy Implant

55, rue Uranus - Z.A.C. Altaïs - 74650 CHAVANOD (France) - Tel.: +33 (0) 450 450 498


@ support.technique@visyimplant.com  victoryimplant.com / visyimplant.com

MADE IN FRANCE

TECHNICAL SHEET Ref. 3.75-MF1

DRILLING PROTOCOL

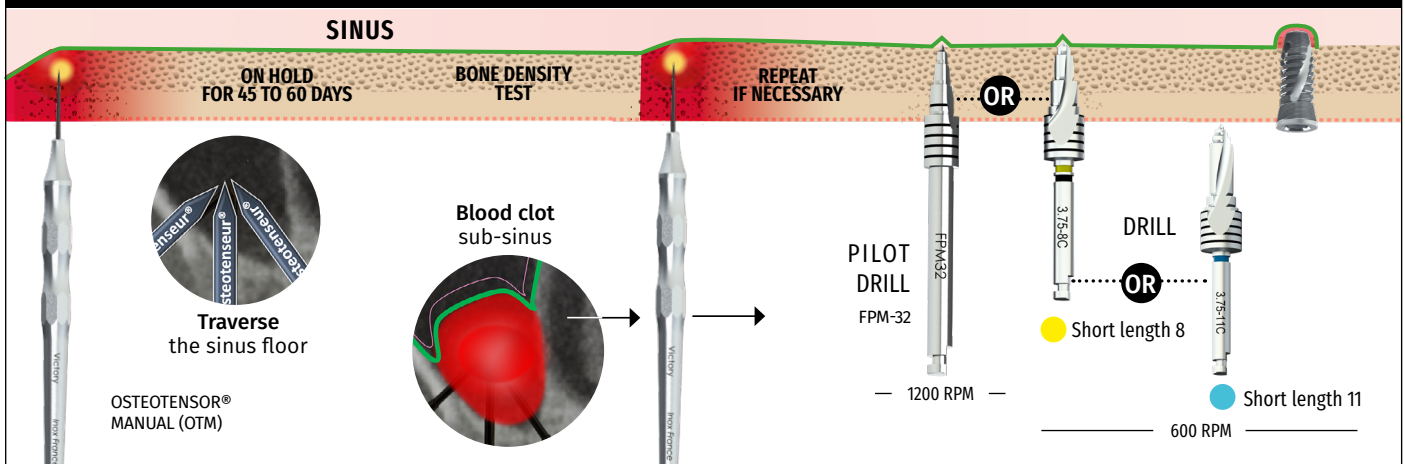
Drilling recommendation: Under continuous irrigation.
Drilling torque recommendation: 5 to 10 N.cm. If the drill bit stalls at 10 N.cm, plan to replace it.

⚠ The drill tip has a length of 0.75 mm. This value must be taken into account during drilling in addition to the implant length. The values D1 to D4 indicate the bone densities for which the drill can be used. Do not use a rotational speed (in revolutions per minute) outside the indicated range.

Implant tightening recommendation: 30 RPM without irrigation.


PLACEMENT PROTOCOL

«No Touch» System
 Breakable implant and cover screw

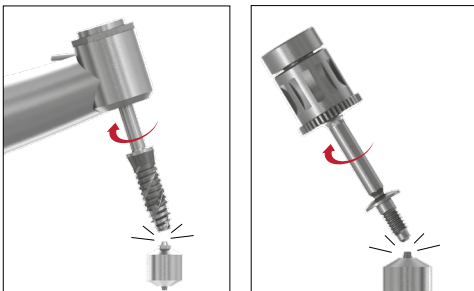
CONTRA-ANGLE

MEC or MEL
 Contra-angle placement

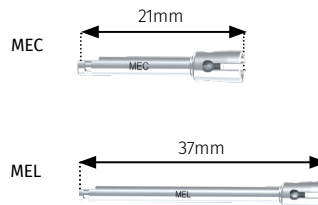
MANUAL

CMC-15 or CMC-30
 Manual placement

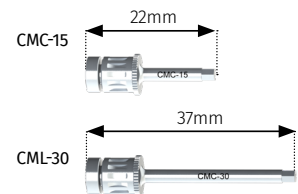TRACEABILITY LABELS

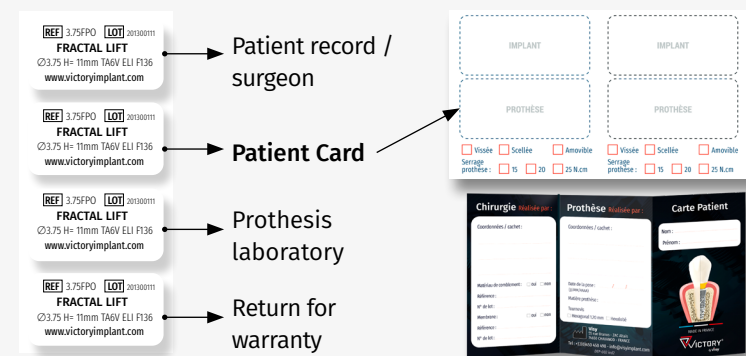CERTIFICATIONS & STANDARDS

Visy is **NF EN ISO 13485 : 2016** certified and holder of **EC certificate** delivered by its notified body for the **93/42/CE** directive. All products are subject to change without notice.


VICTORY WARRANTY


- Victory® products have a 10 years warranty, free from manufacturing defects.
- The buyer remains solely responsible for the use he makes of the products.
- Every component replacements are left to Visy discretion.

RÉFÉRENCES	$\emptyset$ A	$\emptyset$ B	$\emptyset$ C	L
3.3H10-OFX	3.60mm	3.30mm	1.00mm	10mm
3.3H12-OFX	3.60mm			12mm

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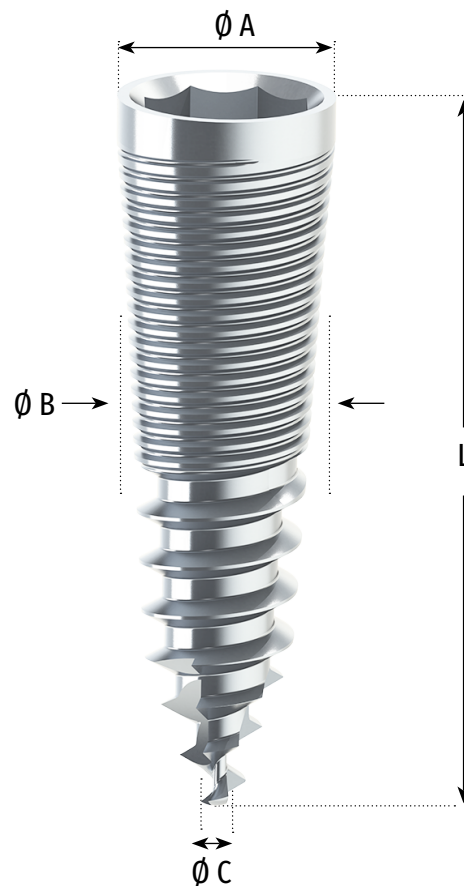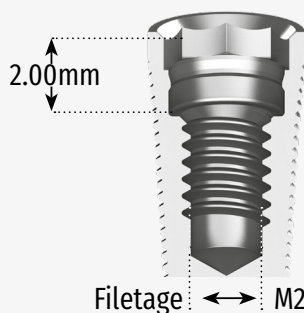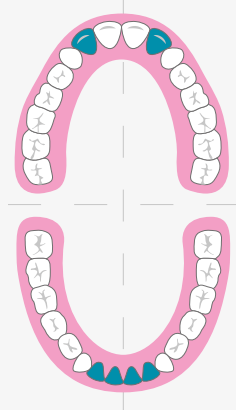
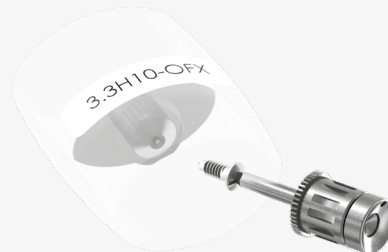
 **MADE IN FRANCE**

LONGUEUR 12


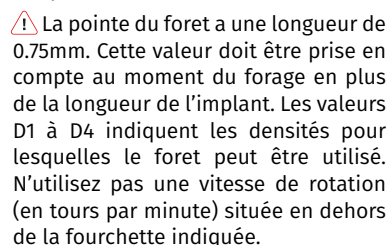
Diagram illustrating the dimensions of the couplers:

- CMC-OI (COURT):** 21mm
- CML-OI (LONG):** 30mm

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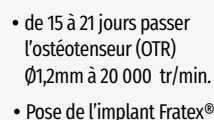Garantie
10 ANS
Victory Implant

Fratex $\emptyset$ 3.30


VICTORY[®]

By Visy


TECHNICAL SHEET Ref. 3.30-OFX-H10 & H12

H10 & H12

REFERENCES	$\emptyset$ A	$\emptyset$ B	$\emptyset$ C	L
3.3H10-OFX	3.60mm	3.30mm	1.00mm	10mm
3.3H12-OFX	3.60mm			12mm

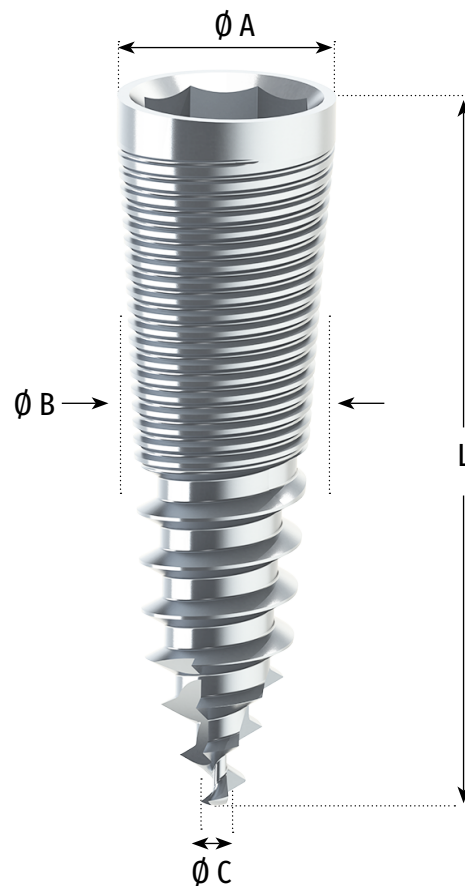INTERNAL OCTAGON PLATFORM

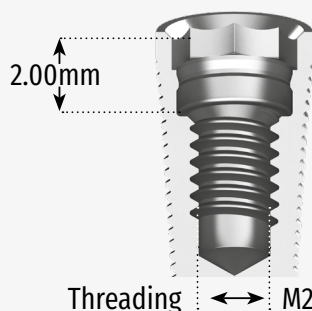- Indexed flat connection by a 2.6mm hexagon
- Allows axis correction up to 90° between two implants

Material


Titanium alloy TA6V ELI F136 - EN ISO 5832-3/ASTM F136

Surface state

VMSS (Victory Smooth Surface)

Roughness: Machined smooth area $0.6 \mu\text{m} < SA < 0.8 \mu\text{m}$


INDICATIONS

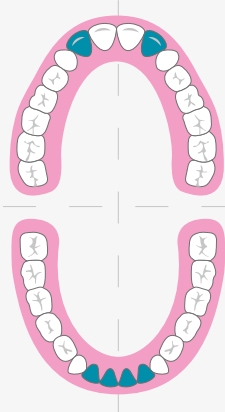- For the replacement of lateral incisors and mandibular incisors
- Prothesis stabilization removable or fixed (unitary or interleaved). The volume and bone quality must be sufficient to support the implants.

SUPPLIED AS STANDARD WITH

- Internal octagon $\emptyset$ 3.60mm breakable cover screw

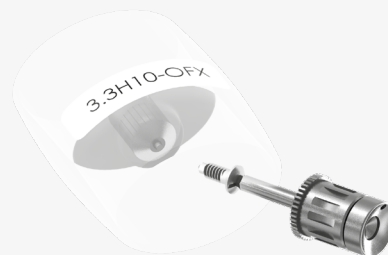Visy Implant

55, rue Uranus - Z.A.C. Altaïs - 74650 CHAVANOD (France) - Tel.: +33 (0) 450 450 498

@ support.technique@visyimplant.com victoryimplant.com / visyimplant.com


victoryimplant.com / visyimplant.com

MADE IN FRANCE

DRILLING PROTOCOL

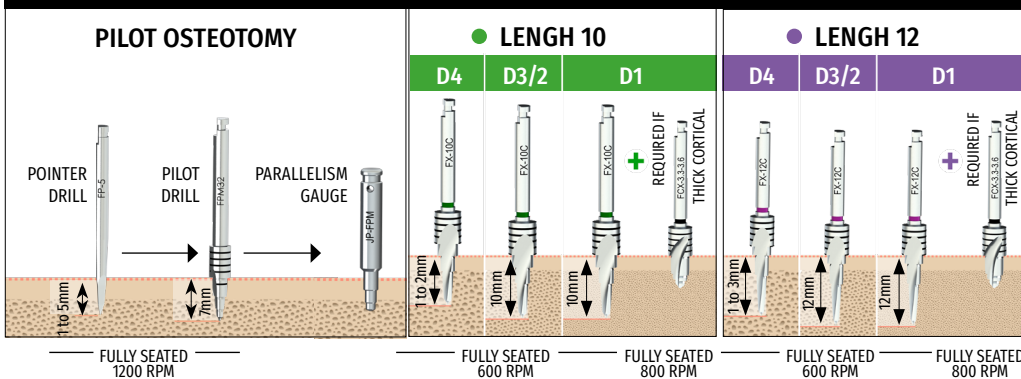Drilling recommendation:

Under continuous irrigation.

Drilling torque recommendation:
5 to 10 N.cm.


If the drill bit stalls at 10 N.cm, plan to replace it.

⚠ The tip of the drill has a length of 0.75 mm. This value must be taken into account during drilling, in addition to the length of the implant. The values D1 to D4 indicate the bone densities for which the drill can be used. Do not use a rotation speed (in revolutions per minute) outside the specified range.

PLACEMENT PROTOCOL

«No Touch» System


Breakable implant and cover screw

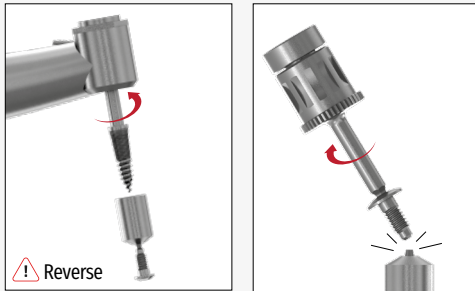CONTRA-ANGLE

MCC-OI or MCL-OI


Contra-angle placement

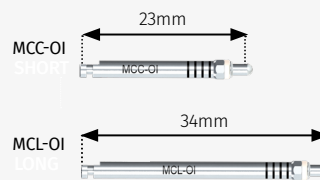MANUAL

CMC-OI or CML-OI

Manual placement

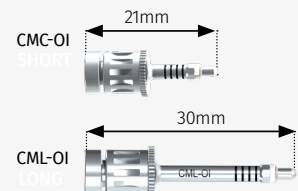Implant drilling recommendation: 30 RPM without irrigation.


1 OSTEOGENIC PREPARATION

- Mark the bone crest
- Use the osteotensor (OTR) Ø1.2 mm at 20,000 RPM

2 OSTEOGENIC ACTIVATION

- 7 to 21 days for transformation of D1 bone into active D2 bone
- Softening of the matrix

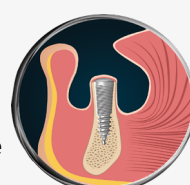
3 EXPANSION OF NARROW RIDGE & IMPLANT PLACEMENT

- Between 15 and 21 days, use the osteotensor (OTR) Ø1.2 mm at 20,000 RPM
- Placement of the Fratex® implant

⚠ Wait 4 months before loading

TRACEABILITY LABELS

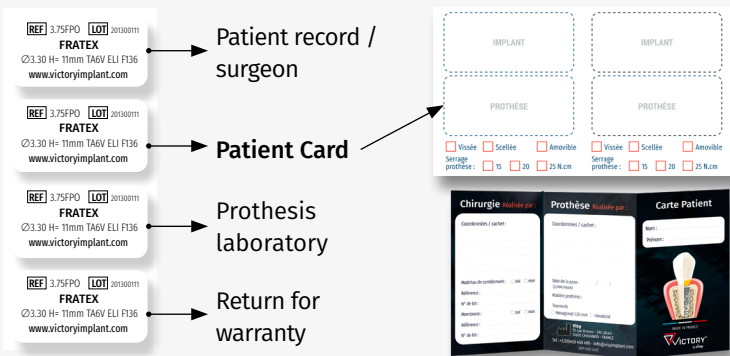CERTIFICATIONS & STANDARDS

Visy is **NF EN ISO 13485 : 2016** certified and holder of **EC certificate** delivered by its notified body for the **93/42/CE** directive. All products are subject to change without notice.

VICTORY WARRANTY


- Victory® products have a 10 years warranty, free from manufacturing defects.
- The buyer remains solely responsible for the use he makes of the products.
- Every component replacements are left to Visy discretion.

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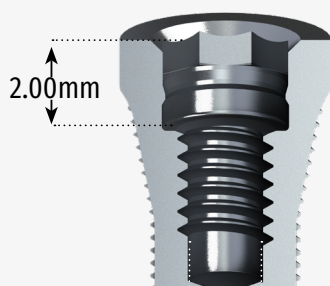
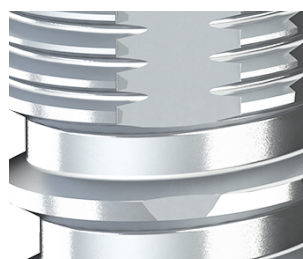
Roughness : Machined smooth area $0.6 \mu\text{m} < SA < 0.8 \mu\text{m}$


INDICATIONS

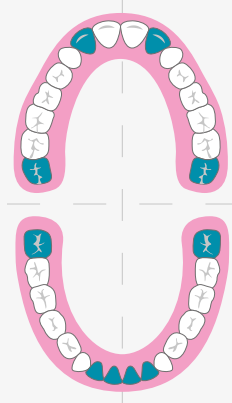- For replacement in the tubero-ptyergoid region, maxillary lateral incisors, and mandibular incisors.
- Prothesis stabilization removable or fixed (unitary or interleaved). The volume and bone quality must be sufficient to support the implants.

SUPPLIED A STARNDARD WITH

- **Internal octagon $\emptyset$ 4.50mm** breakable cover screw

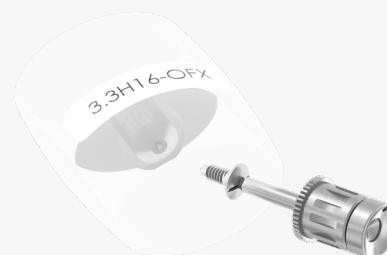Visy Implant

55, rue Uranus - Z.A.C. Altaïs - 74650 CHAVANOD (France) - Tel.: +33 (0) 450 450 498

@ support.technique@visyimplant.com  victoryimplant.com / visyimplant.com

MADE IN FRANCE


Fratex Ø 3.30 H16

TECHNICAL SHEET Ref. 3.3H16-OFX

DRILLING PROTOCOL

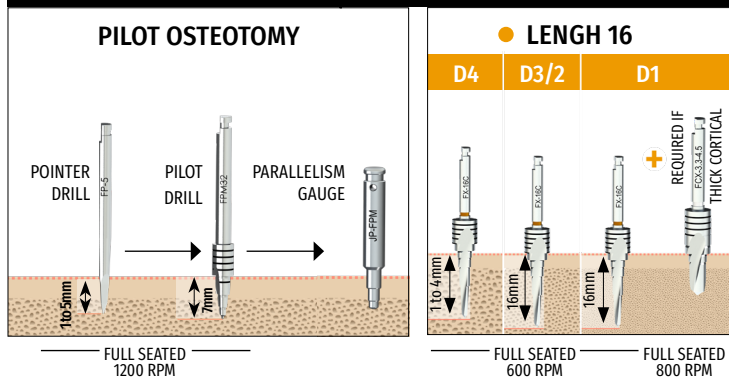
i Drilling protocol: Under continuous irrigation.
Drilling torque recommendation: 5 to 10 N.cm.
If the drill bit stalls at 10 N.cm, plan to replace it.

⚠ The tip of the drill has a length of 0.75 mm. This value must be taken into account during drilling, in addition to the implant length. The values D1 to D4 indicate the bone densities for which the drill can be used. Do not use a rotation speed (in revolutions per minute) outside the indicated range.

Recommended implant screwing protocol: 30 RPM without irrigation.


PLACEMENT PROTOCOL

«No Touch» System
Breakable implant and cover screw

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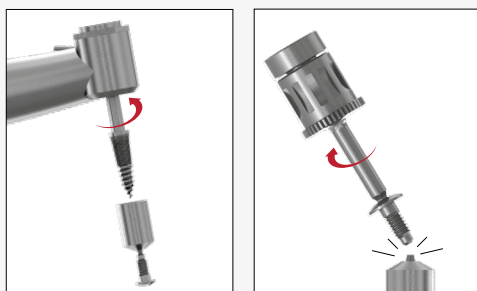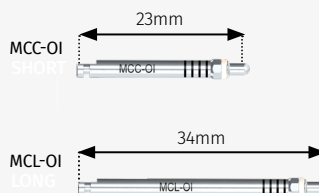
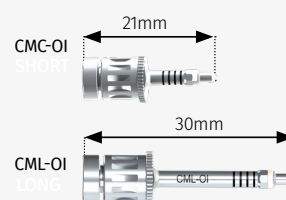
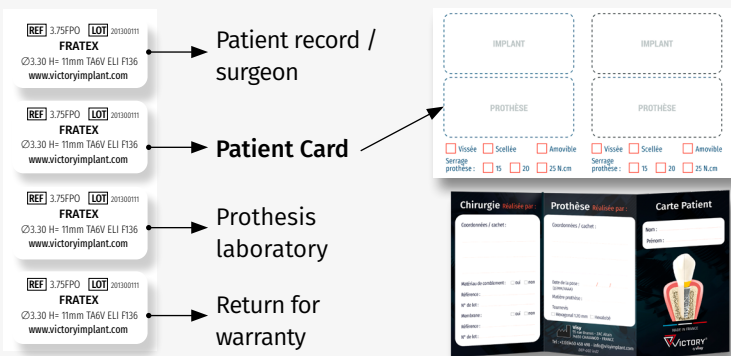
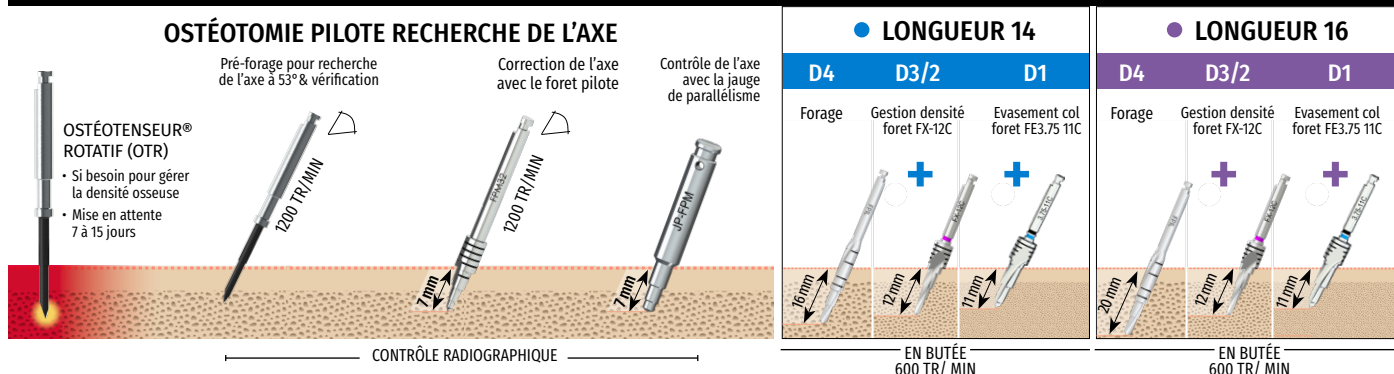
VICTORY
By Visy


TECHNICAL SHEET Ref. 3.30-MFX

REFERENCES	Ø A	Ø B	Ø C	L	L2
3.3H14-MFX	4.50mm	3.30mm	1.00mm	14mm	4mm
3.3H18-MFX				18mm	4mm

MONOBLOC PLATFORM

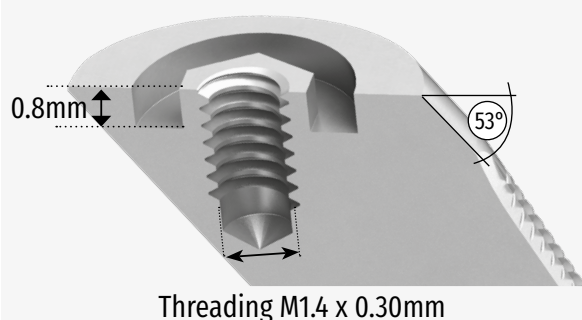- Flat connection Indexed by a 2mm hexagon
- Allows axis correction up to 90° between two implants

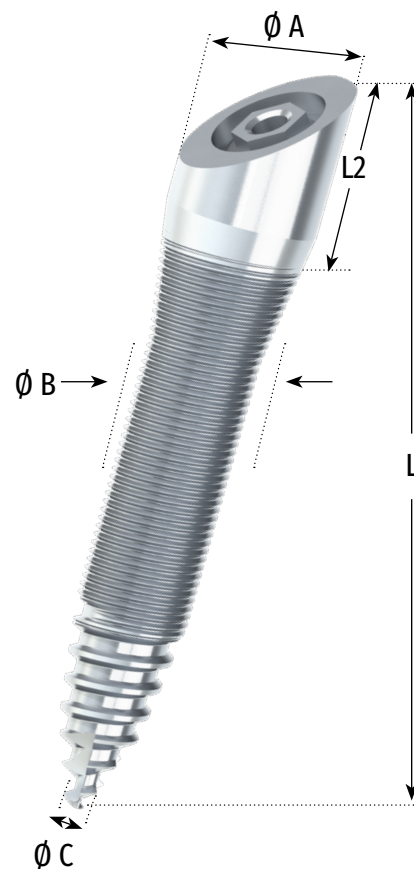Material


Titanium alloy TA6V ELI F136 - EN ISO 5832-3/ASTM F136

Surface state

VMSS (Victory Smooth Surface)


Roughness : Machined smooth area $0.6 \mu\text{m} < SA < 0.8 \mu\text{m}$

INDICATIONS

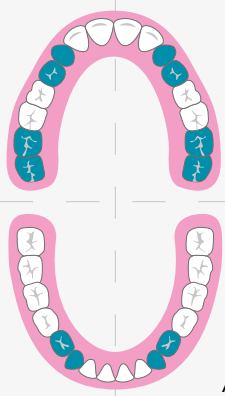
- For the replacement in the sub-sinus and chin region
- Prosthesis stabilization removable or fixed (unitary or interleaved). The volume and bone quality must be sufficient to support the implants.

Adapted to the All-on-Five technique

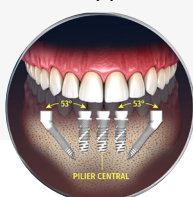SUPPLIED AS STANDARD WITH

- **MONOBLOC Ø 4.50mm** breakable cover screw

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

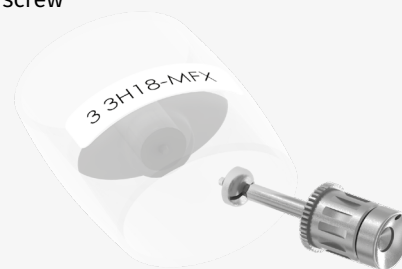
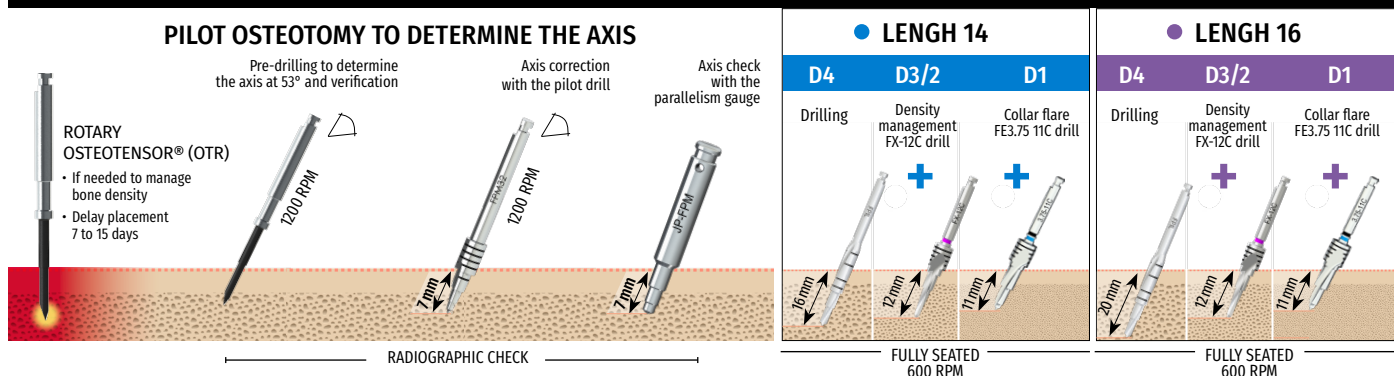
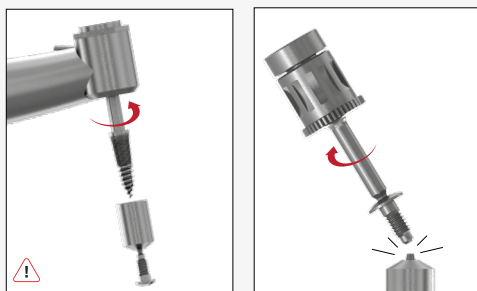
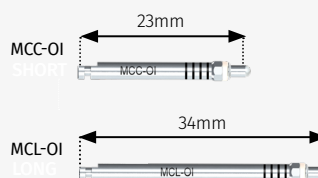
MANUAL

CMC-OI or CML-OI
Manual placement

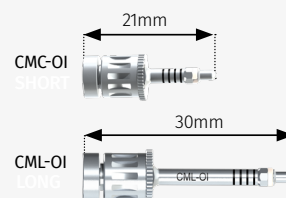Recommended implant insertion speed: 30 rpm without irrigation.

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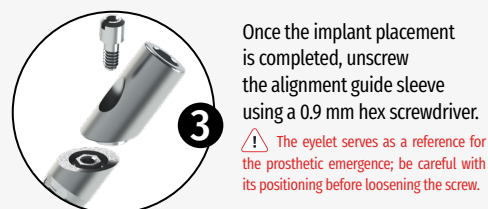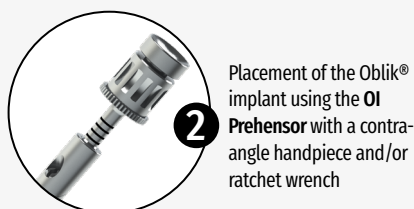
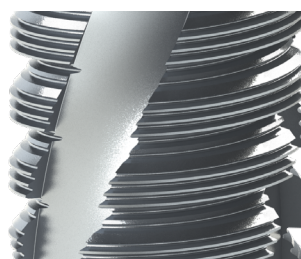
Roughness : Machined smooth area $0.6 \mu\text{m} < SA < 0.8 \mu\text{m}$

INDICATIONS

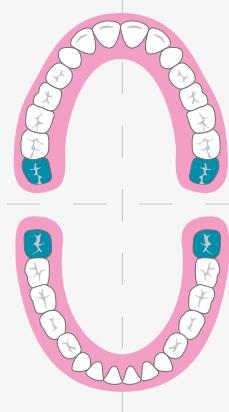- For the replacement in the tubero-ptyergoid region
- Prothesis stabilization removable or fixed (unitary or interleaved). The volume and bone quality must be sufficient to support the implants.

SUPPLIED AS STANDARD WITH

- **MONOBLOC Ø 4.50mm** self breaking cover screw

Visy Implant

55, rue Uranus - Z.A.C. Altaïs - 74650 CHAVANOD (France) - Tel.: +33 (0) 450 450 498

@ support.technique@visyimplant.com victoryimplant.com / visyimplant.com

victoryimplant.com / visyimplant.com

MADE IN FRANCE

PTERYGOID Ø 3.30

TECHNICAL SHEET Ref. 3.30-MF4

DRILLING PROTOCOL

1


Manually insert the Osteotensor® into the tuberosity and verify the axis using a periapical X-ray. The axis should be oriented $\approx 20^\circ/30^\circ$ in the horizontal plane and $\approx 45^\circ$ relative to the maxillary plane.

2

Maintain the axis of the Osteotensor® and gradually guide the pilot drill (Ref. FPIXL) into position (max 600 RPM).

3

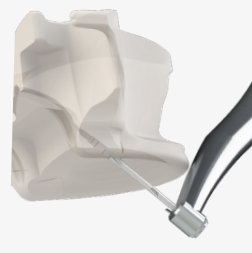
Drilling recommendation: Under continuous irrigation
Recommended drilling torque: 5 to 10 N.cm.
If the drill bit stalls at 10 N-cm, plan to replace it.

PLACEMENT PROTOCOL

«No Touch» System

Breakable implant and cover screw

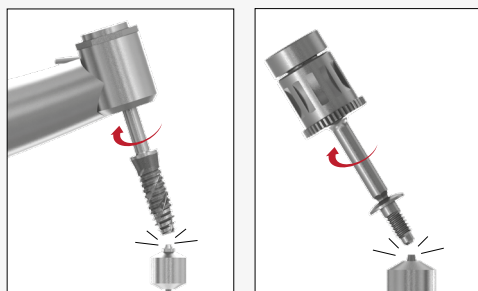4


Implant placement recommendation:
Speed : 30 RPM without irrigation
Torque : 35 to 65 N.cm

CONTRA-ANGLE

MEC or MEL


Contra-angle placement

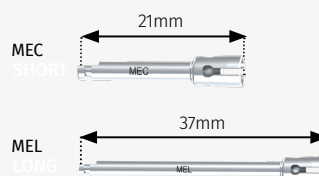MANUAL


CMC-15 or CMC-30

Manual placement

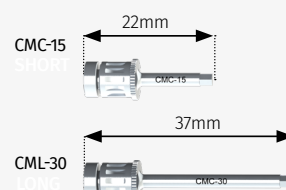5

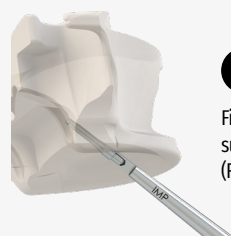Firmly impact the implant with a surgical mallet and the impactor (Ref. IMT).

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