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FICHES TECHNIQUES - TECHNICAL SHEET

FICHES TECHNIQUES

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Fractal Plus	Ø 3.75 FPO	2	
	Ø 4.75 FPO	6	
Fractal Lift	Ø 3.75 OF1	10	
	Ø 4.75 OF1	14	
	Ø 3.75 MF1	18	
Fratex	Ø 3.30 H10 & H12 OFX	22	
	Ø 3.30 H16 OFX	26	
Fratex Oblik	Ø 3.30 MFX	30	
Pterygoïdien	Ø 4.75 MF4	34	

TECHNICAL SHEET

DOWNLOAD ALL 

Fractal Plus	Ø 3.75 FPO	4	
	Ø 4.75 FPO	8	
Fractal Lift	Ø 3.75 OF1	12	
	Ø 4.75 OF1	16	
	Ø 3.75 MF1	20	
Fratex	Ø 3.30 H10 & H12 OFX	24	
	Ø 3.30 H16 OFX	28	
Fratex Oblik	Ø 3.30 MFX	32	
Pterygoïdien	Ø 4.75 MF4	36	

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

INDICATIONS

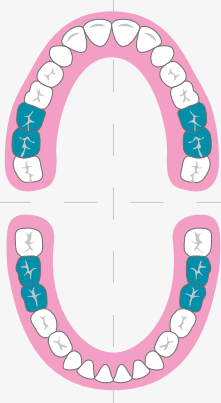- For the mollars replacement
- Prosthesis stabilization removable or fixed (unitary or interleaved). The volume and bone quality must be sufficient to support the implants.

SUPPLIED AS STANDARDS WITH

- Internal octagon $\emptyset$ 4.50mm breakable cover screw

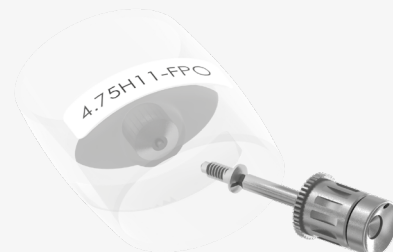Visy Implant


55, rue Uranus - Z.A.C. Altaïs - 74650 CHAVANOD (France) - Tel.: +33 (0) 450 450 498

@ support.technique@visyimplant.com victoryimplant.com / visyimplant.com

MADE IN FRANCE

TECHNICAL SHEET Ref. 4.75-FPO

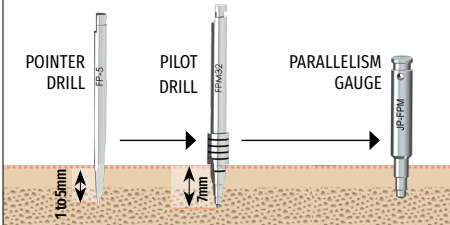DRILLING PROTOCOL		● LENGTH 8			● LENGTH 11		
PILOT OSTEOTOMY		D4	D3/2	D1	D4	D3/2	D1
				 + REQUIRED IF FCF-4.75 THICK CORTICAL			 + REQUIRED IF FCF-4.75 THICK CORTICAL
IN STOP. 1200 RPM		IN STOP. 600 RPM			IN STOP. 600 RPM		
		IN STOP. 800 RPM			IN STOP. 800 TR/ MIN		

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

INDICATIONS

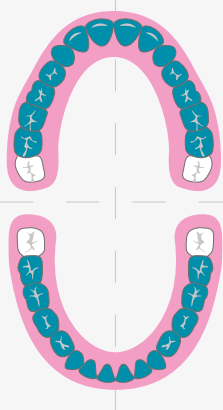- For the replacement in all sectors
- **Prosthesis stabilization** removable or fixed (unitary or interleaved). The volume and bone quality must be sufficient to support the implants.

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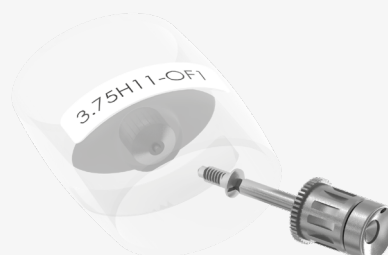
MADE IN FRANCE

DRILLING PROTOCOL

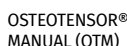Drilling torque recommendation: 5 to 10 N.cm. If the drill bit stalls at 10 N-cm, plan to replace it.

The values D1 to D4 indicate the bone densities for which the drill can be used. Do not use a rotational speed (in revolutions per minute) outside the indicated range.

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

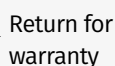
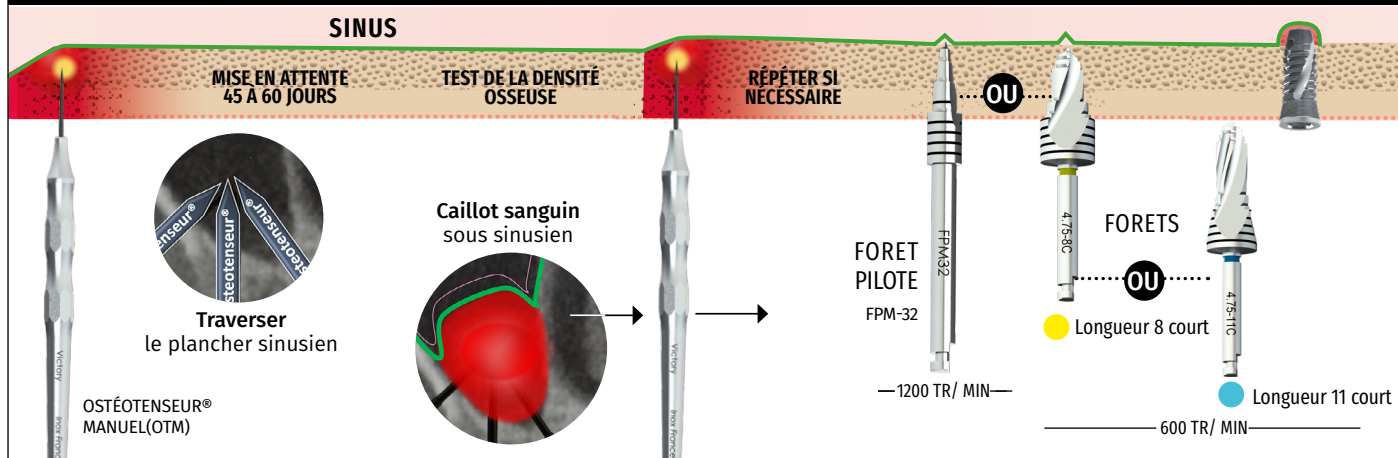
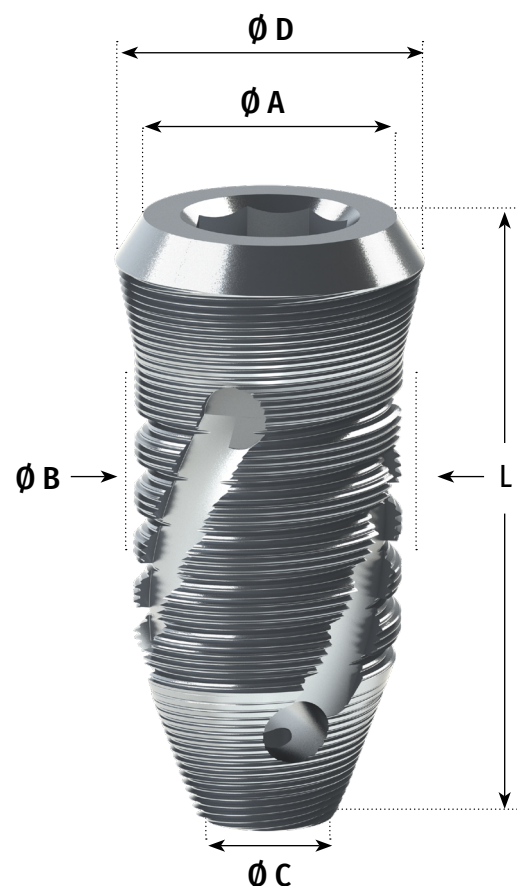
INTERNAL OCTAGON PLATFORM

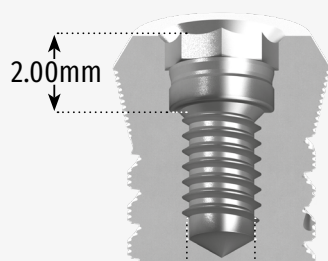Threading $\longleftrightarrow$ M2 x 0.40mm

- Flat indexed connection by a 2.6mm octagon
- Retention groove for centering prosthetic components

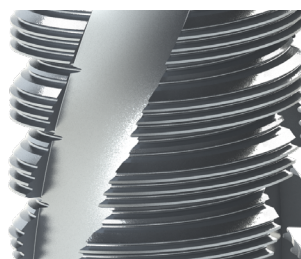Material


Titanium alloy TA6V ELI F136 - EN ISO 5832-3/ASTM F136

Surface state

VMSS (Victory Smooth Surface)

Roughness : Machined smooth area $0.6 \mu\text{m} < SA < 0.8 \mu\text{m}$

INDICATIONS

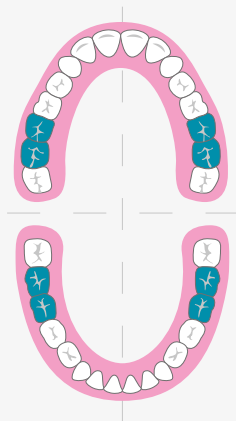- For the mollars replacement
- **Prosthesis stabilization** removable or fixed (unitary or interleaved). The volume and bone quality must be sufficient to support the implants.

SUPPLIED AS STANDARDS WITH

- **Internal octagon $\varnothing$ 4.50mm** breakable cover screw

Visy Implant

55, rue Uranus - Z.A.C. Altaïs - 74650 CHAVANOD (France) - Tel.: +33 (0) 450 450 498


@ support.technique@visyimplant.com  victoryimplant.com / visyimplant.com

MADE IN FRANCE

TECHNICAL SHEET Ref. 4.75-OF1

DRILLING PROTOCOL

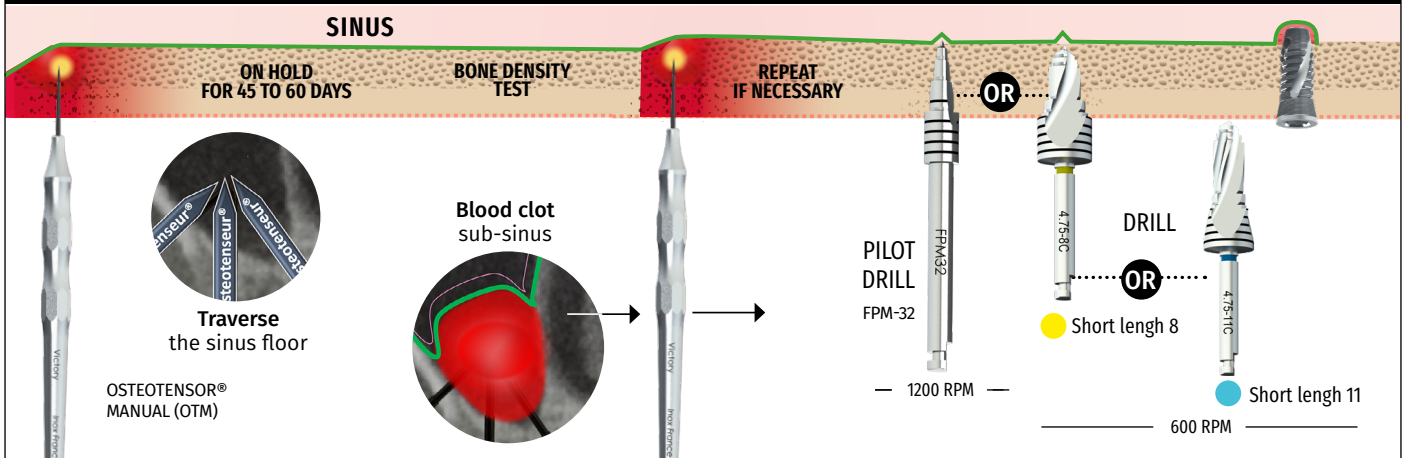
Drilling recommendation: Under continuous irrigation.
Drilling torque recommendation: 5 to 10 N.cm. If the drill bit stalls at 10 N.cm, plan to replace it.

The drill tip has a length of 0.75 mm. This value must be taken into account during drilling in addition to the implant length. The values D1 to D4 indicate the bone densities for which the drill can be used. Do not use a rotational speed (in revolutions per minute) outside the indicated range.

Implant tightening recommendation: 30 RPM without irrigation.


PLACEMENT PROTOCOL

«No Touch» System
 Breakable implant and cover screw

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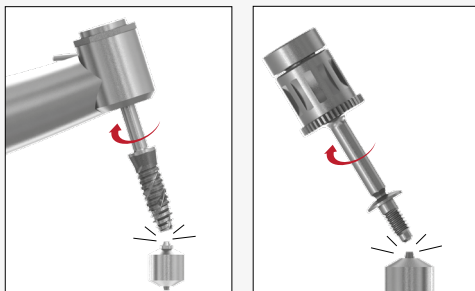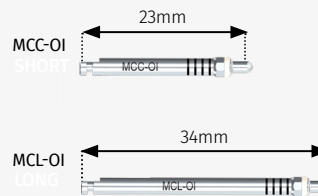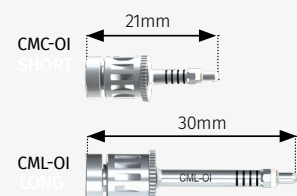
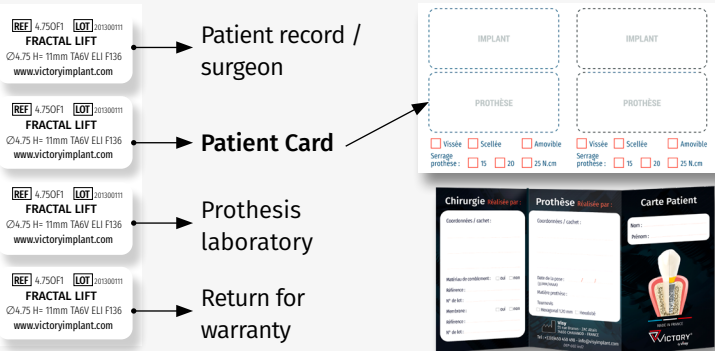
VICTORY WARRANTY

- Victory[®] products have a 10 years warranty, free from manufacturing defects.
- The buyer remains solely responsible for the use he makes of the products.
- Every component replacements are left to Visy discretion.

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


MONOBLOC PLATFORM

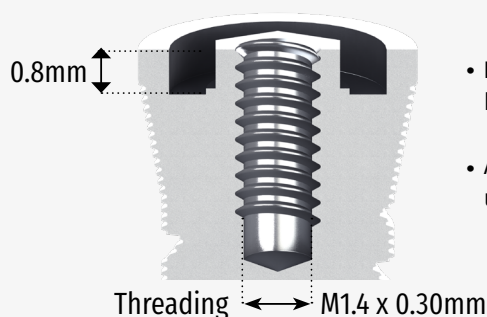- **Indexed flat connection** by a 2mm hexagon
- Allows axis correction up to 90° between two implants

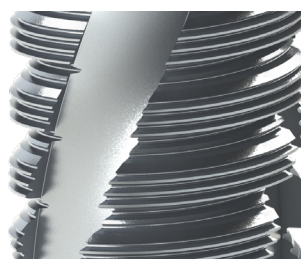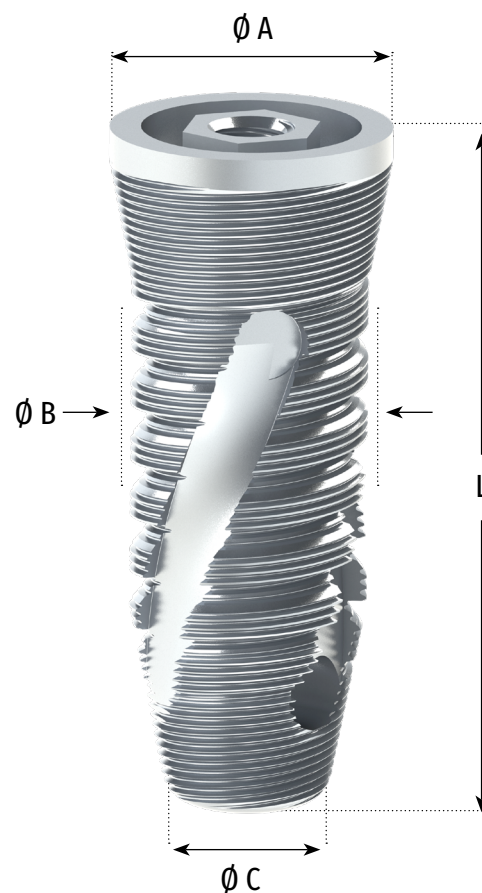Material


Titanium alloy TA6V ELI F136 - EN ISO 5832-3/ASTM F136

Surface state

VMSS (Victory Smooth Surface)

Roughness : Machined smooth area $0.6 \mu\text{m} < SA < 0.8 \mu\text{m}$

INDICATIONS

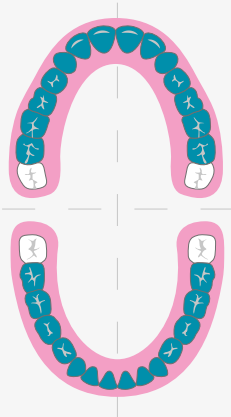- For all sectors replacement
- **Prosthesis stabilization** removable or fixed (unitary or interleaved). The volume and bone quality must be sufficient to support the implants.

SUPPLIED AS STANDARDS WITH

- **MONOBLOC $\varnothing$ 4.50mm** self breaking cover screw

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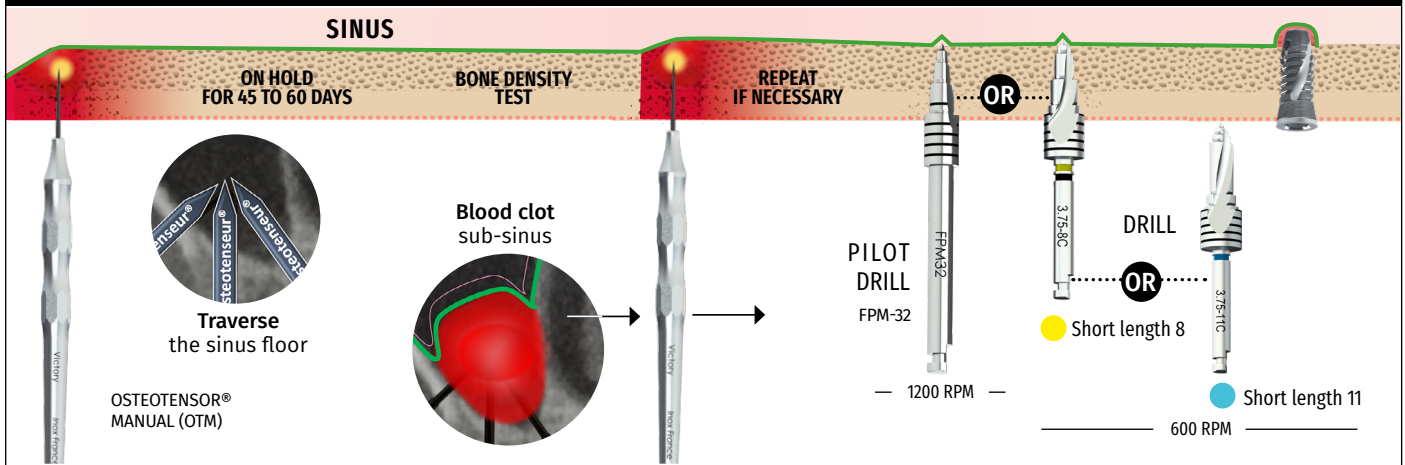
⚠ The drill tip has a length of 0.75 mm. This value must be taken into account during drilling in addition to the implant length. The values D1 to D4 indicate the bone densities for which the drill can be used. Do not use a rotational speed (in revolutions per minute) outside the indicated range.

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
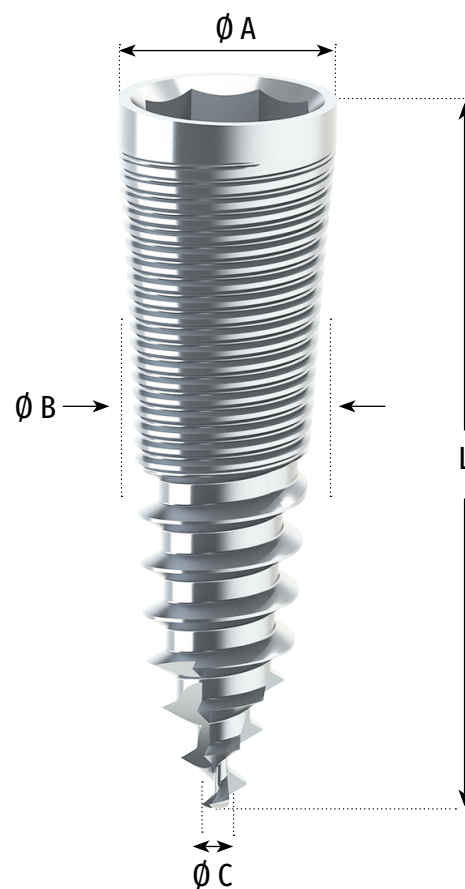
INTERNAL OCTAGON PLATFORM

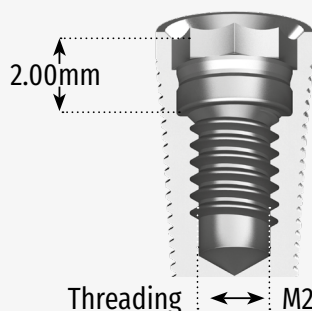- Indexed flat connection by a 2.6mm hexagon
- Allows axis correction up to 90° between two implants

Material


Titanium alloy TA6V ELI F136 - EN ISO 5832-3/ASTM F136

Surface state

VMSS (Victory Smooth Surface)

Roughness: Machined smooth area $0.6 \mu\text{m} < SA < 0.8 \mu\text{m}$


INDICATIONS

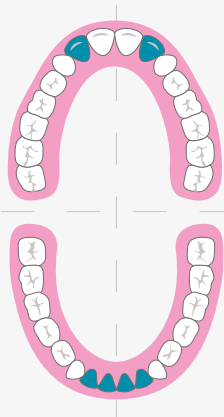- For the replacement of lateral incisors and mandibular incisors
- Prosthesis stabilization removable or fixed (unitary or interleaved). The volume and bone quality must be sufficient to support the implants.

SUPPLIED AS STANDARD WITH

- Internal octagon $\emptyset$ 3.60mm breakable cover screw

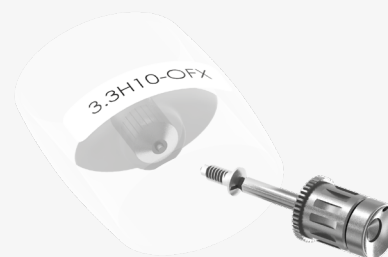Visy Implant

55, rue Uranus - Z.A.C. Altaïs - 74650 CHAVANOD (France) - Tel.: +33 (0) 450 450 498

@ support.technique@visyimplant.com victoryimplant.com / visyimplant.com


victoryimplant.com / visyimplant.com

MADE IN FRANCE

DRILLING PROTOCOL

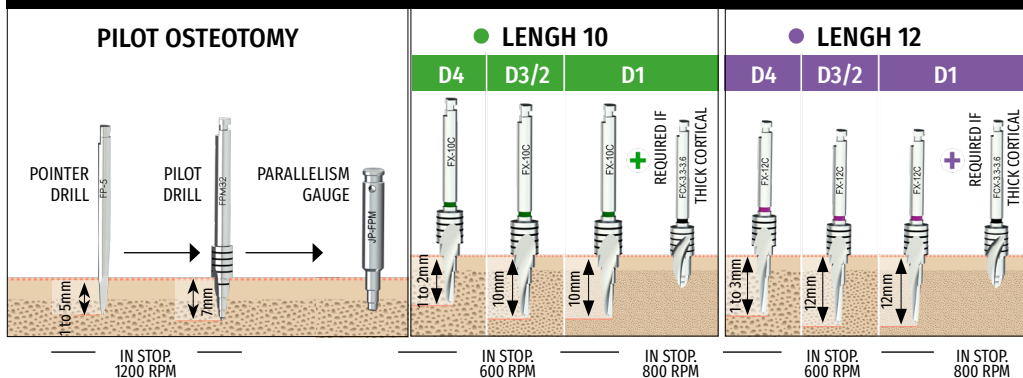Drilling recommendation:

Under continuous irrigation.

Drilling torque recommendation:

5 to 10 N.cm.


If the drill bit stalls at 10 N.cm, plan to replace it.

⚠ The tip of the drill has a length of 0.75 mm. This value must be taken into account during drilling, in addition to the length of the implant. The values D1 to D4 indicate the bone densities for which the drill can be used. Do not use a rotation speed (in revolutions per minute) outside the specified range.

PLACEMENT PROTOCOL

«No Touch» System


Breakable implant and cover screw

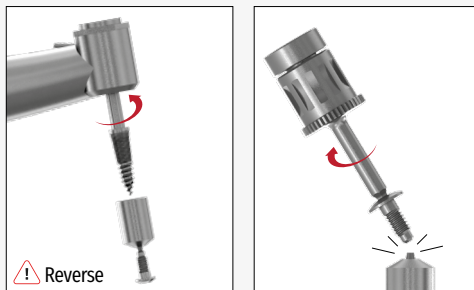CONTRA-ANGLE

MCC-OI or MCL-OI


Contra-angle placement



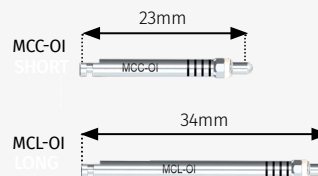
MANUAL

CMC-OI or CML-OI


Manual placement

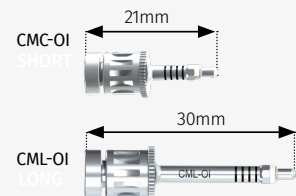Implant drilling recommendation: 30 RPM without irrigation.

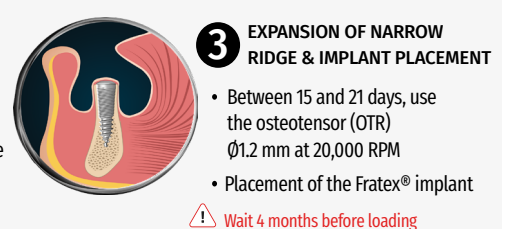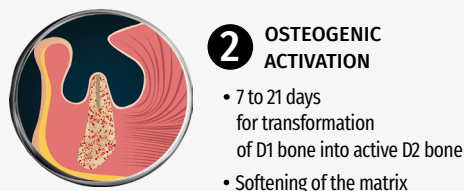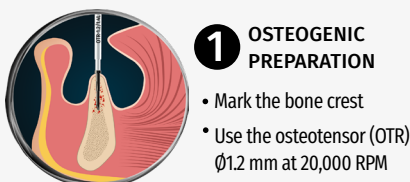TRACEABILITY LABELS

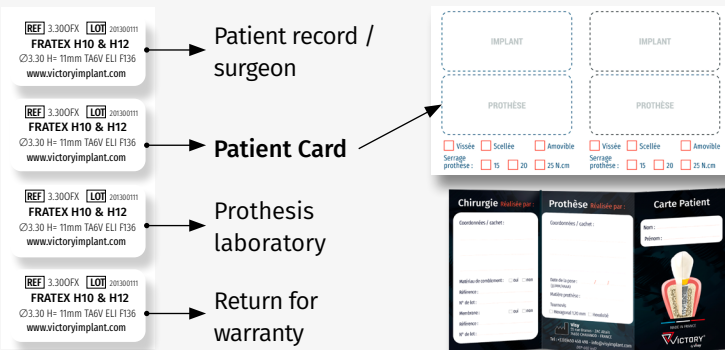CERTIFICATIONS & STANDARDS

Visy is **NF EN ISO 13485 : 2016** certified and holder of **EC certificate** delivered by its notified body for the **93/42/CE** directive. All products are subject to change without notice.


VICTORY WARRANTY

- Victory® products have a 10 years warranty, free from manufacturing defects.
- The buyer remains solely responsible for the use he makes of the products.
- Every component replacements are left to Visy discretion.

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
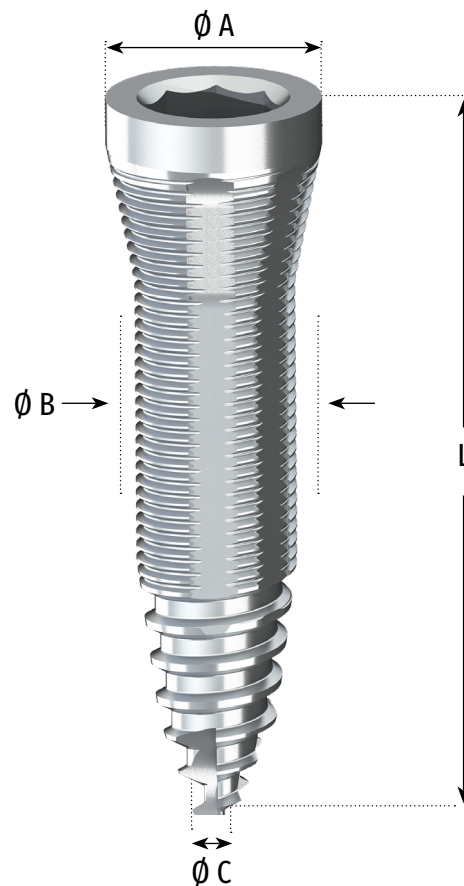
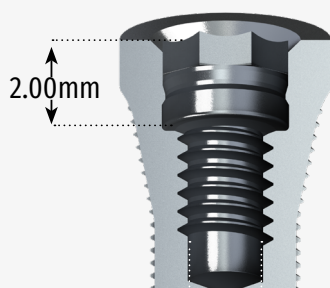
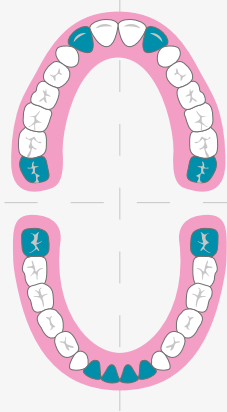
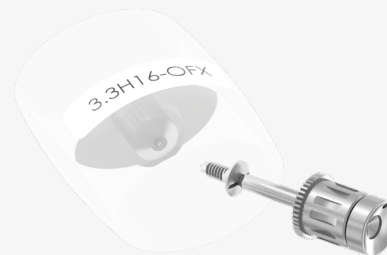
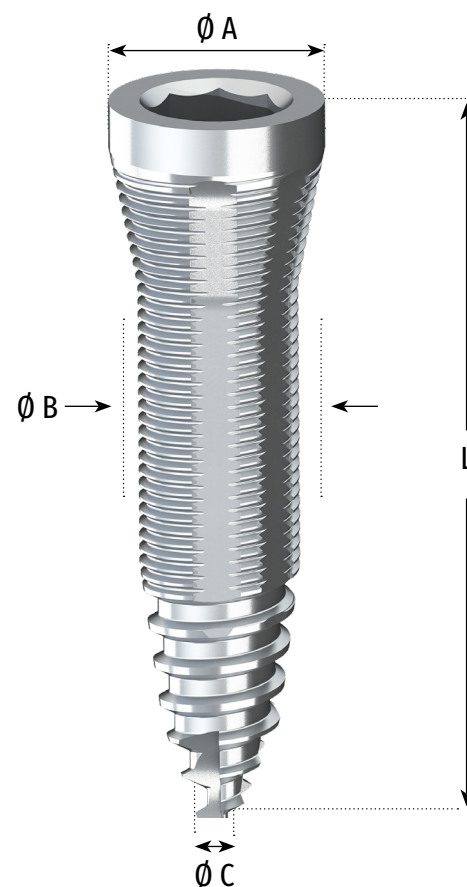
INTERNAL OCTAGON PLATFORM

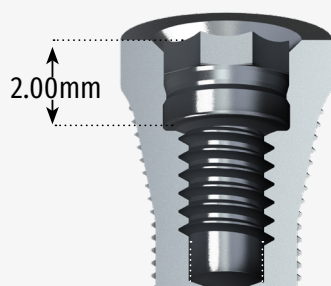Threading $\longleftrightarrow$ M2 x 0.40mm

- Indexed flat connection by a 2.6mm hexagon
- Allows axis correction up to 90° between two implants

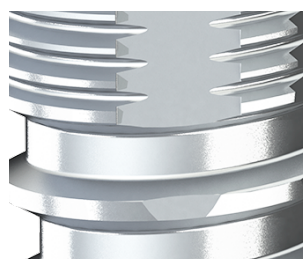Material


Titanium alloy TA6V ELI F136 - EN ISO 5832-3/ASTM F136

Surface state

VMSS (Victory Smooth Surface)

Roughness : Machined smooth area $0.6 \mu\text{m} < SA < 0.8 \mu\text{m}$


INDICATIONS

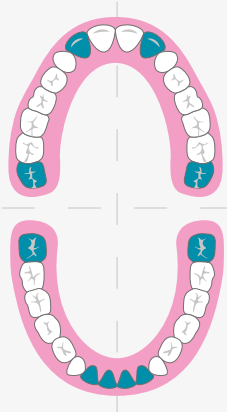- For replacement in the tubero-ptyergoid region, maxillary lateral incisors, and mandibular incisors.
- Prosthesis stabilization removable or fixed (unitary or interleaved). The volume and bone quality must be sufficient to support the implants.

SUPPLIED A STARNDARD WITH

- **Internal octagon $\emptyset$ 4.50mm** breakable cover screw

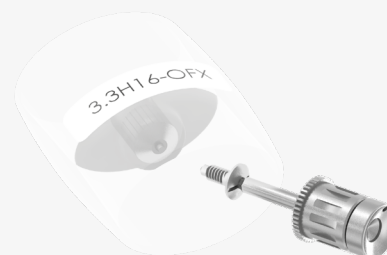Visy Implant

55, rue Uranus - Z.A.C. Altaïs - 74650 CHAVANOD (France) - Tel.: +33 (0) 450 450 498


@ support.technique@visyimplant.com  victoryimplant.com / visyimplant.com

MADE IN FRANCE

TECHNICAL SHEET **Ref. 3.3H16-OFX**

DRILLING PROTOCOL

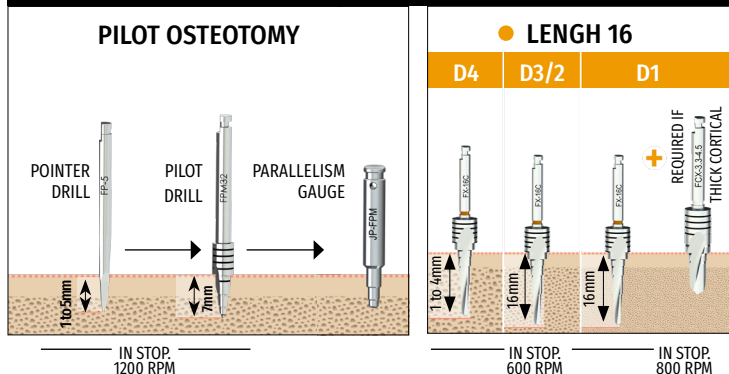i **Drilling protocol:** Under continuous irrigation.
Drilling torque recommendation: 5 to 10 N.cm.
 If the drill bit stalls at 10 N.cm, plan to replace it.


⚠ The tip of the drill has a length of 0.75 mm. This value must be taken into account during drilling, in addition to the implant length. The values D1 to D4 indicate the bone densities for which the drill can be used. Do not use a rotation speed (in revolutions per minute) outside the indicated range.

Recommended implant screwing protocol: 30 RPM without irrigation.

PLACEMENT PROTOCOL

«No Touch» System


Breakable implant and cover screw

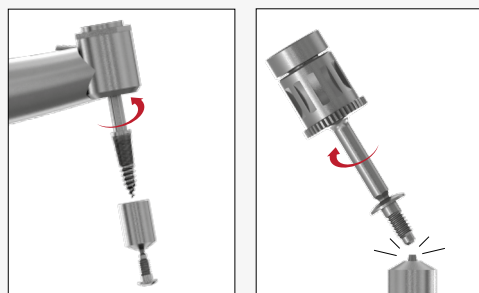CONTRA-ANGLE

MCC-OI or MCL-OI


Contra-angle placement

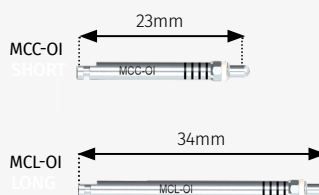
MANUAL

CMC-OI or CML-OI

Manual placement

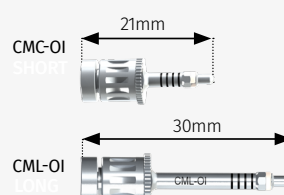TRACEABILITY LABELS

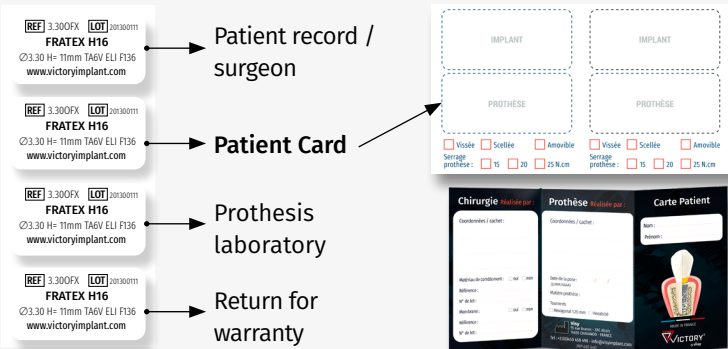CERTIFICATIONS & STANDARDS

Visy is **NF EN ISO 13485 : 2016** certified and holder of **EC certificate** delivered by its notified body for the **93/42/CE** directive. All products are subject to change without notice.


VICTORY WARRANTY

- Victory® products have a 10 years warranty, free from manufacturing defects.
- The buyer remains solely responsible for the use he makes of the products.
- Every composant replacements are left to Visy discretion.

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
RÉFÉRENCES	Ø A	Ø B	Ø C	L	L2
3.3H14-MFX	4.50mm	3.30mm	1.00mm	14mm	4mm
3.3H18-MFX				18mm	4mm

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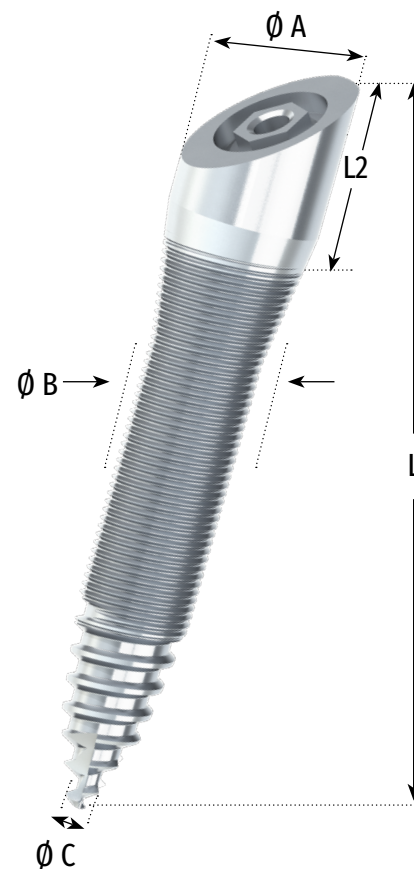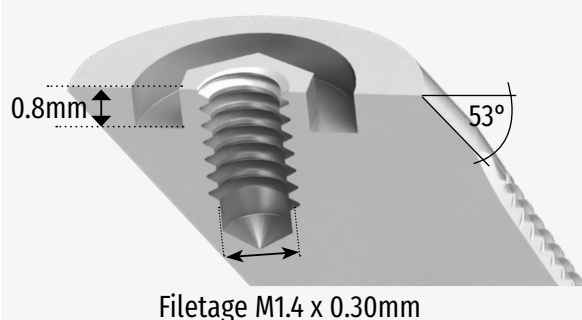
Fratex Oblik Ø 3.30


TECHNICAL SHEET Ref. 3.30-MFX

REFERENCES	Ø A	Ø B	Ø C	L	L2
3.3H14-MFX	4.50mm	3.30mm	1.00mm	14mm	4mm
3.3H18-MFX				18mm	4mm

MONOBLOC PLATFORM

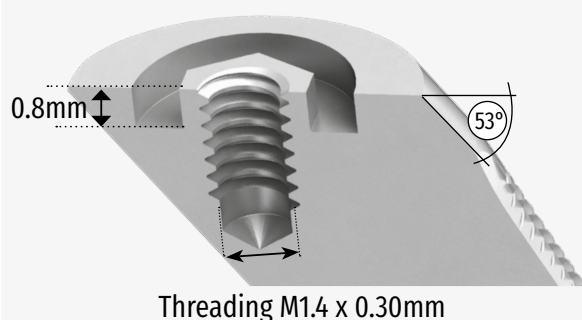- Flat connection Indexed by a 2mm hexagon
- Allows axis correction up to 90° between two implants

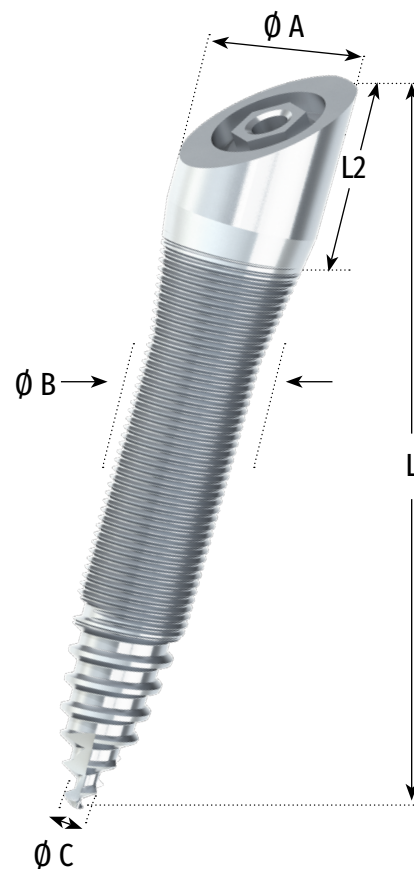Material


Titanium alloy TA6V ELI F136 - EN ISO 5832-3/ASTM F136

Surface state

VMSS (Victory Smooth Surface)


Roughness : Machined smooth area $0.6 \mu\text{m} < SA < 0.8 \mu\text{m}$

INDICATIONS

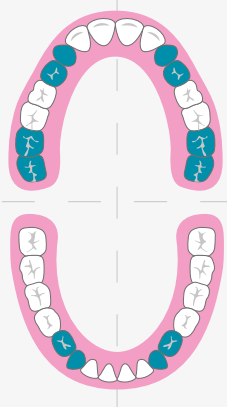
- For the replacement in the sub-sinus and chin region
- Prosthesis stabilization removable or fixed (unitary or interleaved). The volume and bone quality must be sufficient to support the implants.

Adapted to the All-on-Five technique

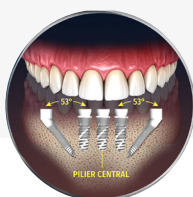SUPPLIED AS STANDARD WITH

- **MONOBLOC Ø 4.50mm** breakable cover screw

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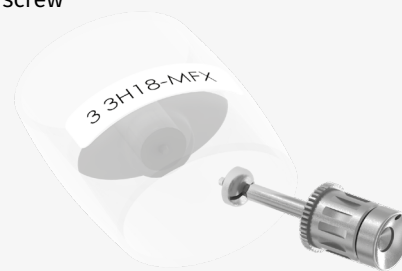
TECHNICAL SHEET **Ref. 3.30-MFX**

DRILLING PROTOCOL


PILOT OSTEOTOMY TO DETERMINE THE AXIS

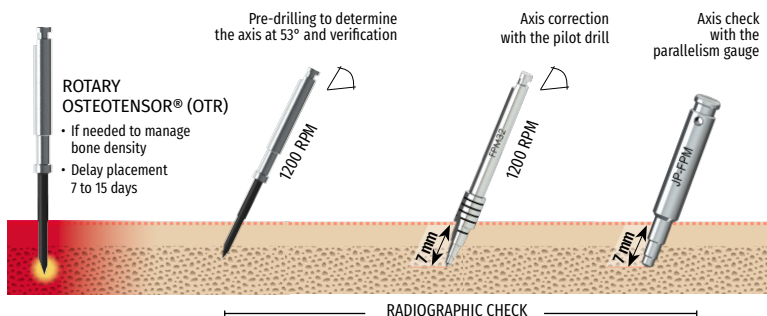● LENGH 14			● LENGH 16		
D4	D3/2	D1	D4	D3/2	D1
Drilling	Density management FX-12C drill	Collar flare FE3.75 11C drill	Drilling	Density management FX-12C drill	Collar flare FE3.75 11C drill
16 mm	12 mm	11 mm	20 mm	12 mm	11 mm
IN STOP. 600 RPM			IN STOP. 600 RPM		


PLACEMENT PROTOCOL

«No Touch» System Breakable implant and cover screw

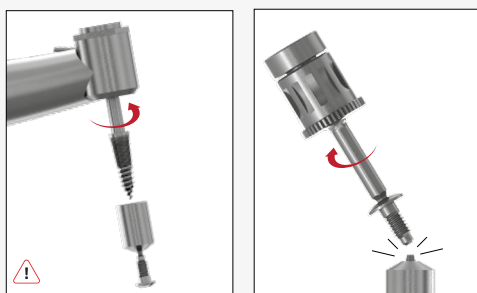
CONTRA-ANGLE

MCC-OI or MCL-OI Contra-angle placement

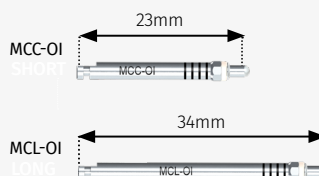MANUAL

CMC-OI or CML-OI Manual placement

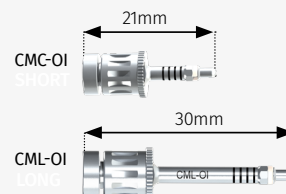Recommended implant insertion speed: 30 rpm without irrigation.

1 Alignment correction guide sleeve, mounted on the Oblik® implant

2 Placement of the Oblik® implant using the OI Prehensor with a contra-angle handpiece and/or ratchet wrench

Once the implant placement is completed, unscrew the alignment guide sleeve using a 0.9 mm hex screwdriver.

⚠ The eyelet serves as a reference for the prosthetic emergence; be careful with its positioning before loosening the screw.

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

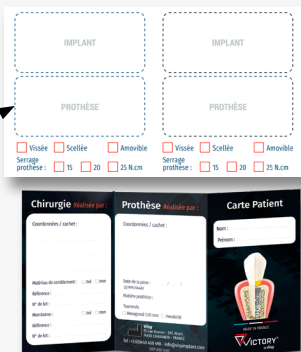
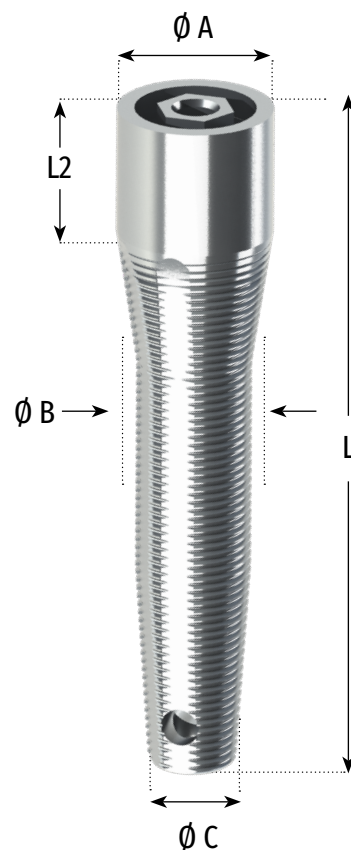
MONOBLOC PLATFORM

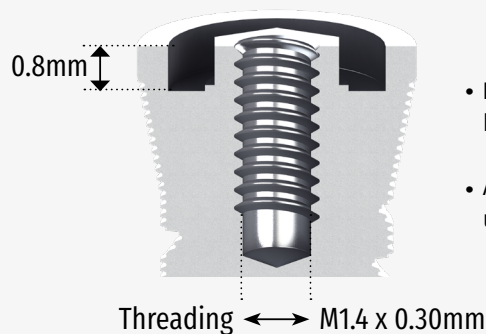- Indexed flat connection by a 2mm hexagon
- Allows axis correction up to 90° between two implants

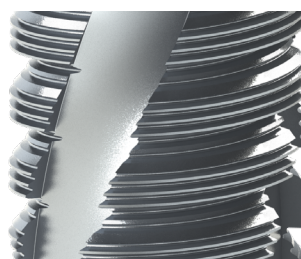Material


Titanium alloy TA6V ELI F136 - EN ISO 5832-3/ASTM F136

Surface state

VMSS (Victory Smooth Surface)

Roughness : Machined smooth area $0.6 \mu\text{m} < SA < 0.8 \mu\text{m}$


INDICATIONS

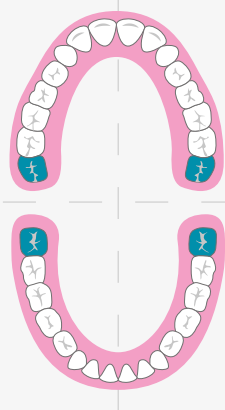- For the replacement in the tubero-ptyergoid region
- Prothesis stabilization removable or fixed (unitary or interleaved). The volume and bone quality must be sufficient to support the implants.

SUPPLIED AS STANDARD WITH

- **MONOBLOC Ø 4.50mm** self breaking cover screw

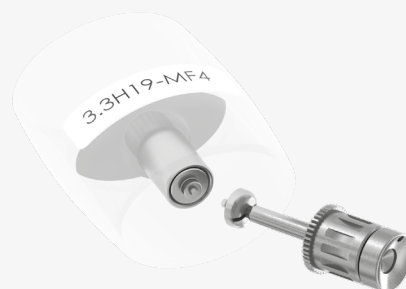Visy Implant

55, rue Uranus - Z.A.C. Altaïs - 74650 CHAVANOD (France) - Tel.: +33 (0) 450 450 498

@ support.technique@visyimplant.com victoryimplant.com / visyimplant.com

victoryimplant.com / visyimplant.com

MADE IN FRANCE

PTERYGOID Ø 3.30

TECHNICAL SHEET Ref. 3.30-MF4

DRILLING PROTOCOL

1


Manually insert the Osteotensor® into the tuberosity and verify the axis using a periapical X-ray. The axis should be oriented $\approx 20^\circ/30^\circ$ in the horizontal plane and $\approx 45^\circ$ relative to the maxillary plane.

2

Maintain the axis of the Osteotensor® and gradually guide the pilot drill (Ref. FPIXL) into position (max 600 RPM).

3

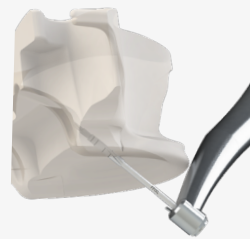
Drilling recommendation: Under continuous irrigation
Recommended drilling torque: 5 to 10 N.cm.
If the drill bit stalls at 10 N-cm, plan to replace it.

PLACEMENT PROTOCOL

«No Touch» System


Breakable implant and cover screw

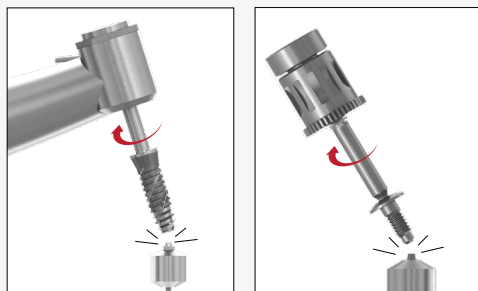CONTRA-ANGLE

MEC or MEL


Contra-angle placement

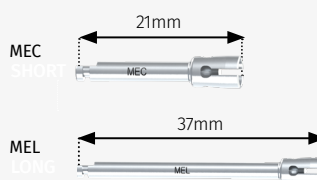MANUAL

CMC-15 or CMC-30

Manual placement

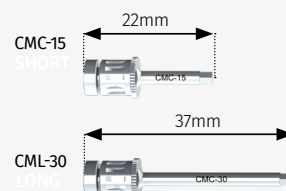
4

Implant placement recommendation:
Speed : 30 RPM without irrigation
Torque : 35 to 65 N.cm

5

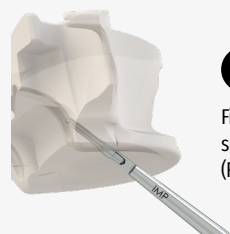Firmly impact the implant with a surgical mallet and the impactor (Ref. IMT).

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